

BEAUTY DEVICES

BEAUTY REJUVENATION DEVICE

User Manual ▶


BEAUTY REJUVENATION DEVICE

D020

User Manual ▶

Please read the instruction carefully before use

PRODUCT STRUCTURE ▶

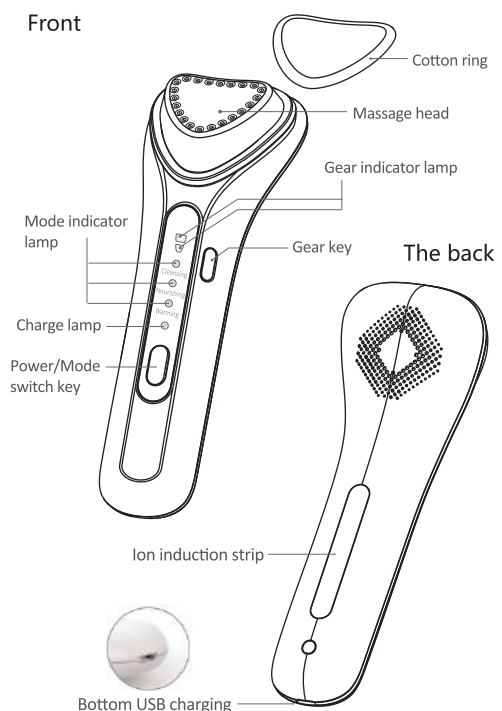

P1

FEATURES ▶

This product is a facial beauty instrument capable of deep cleaning, deep nourishing, warm and moisturizing, and high frequency vibration of positive ion and negative ion.

- Positive and negative ions: by using the principle of positive and negative charge ion heterosexual absorption, deep pore fouling is derived by biological current and thoroughly cleaned. Beauty ingredients are permeated into the muscle floor by ion introduction.
- Vibration: mechanical vibration, voltage stimulation, make the human skin tighten, increase its own elasticity, enhance the metabolism of skin cells.
- Warm heat: warm massage promotes absorption, makes skin texture extend and skin becomes fuller.
- Clean skin pattern: clean the pores deep garbage and lay a good foundation for subsequent skin care.
- Nourishing and importing mode: full introduction of skin care products, enhance cell activity, desalinate fine lines, nourish skin.
- Warm heat moisturizing mode: warm heat massage to promote blood circulation, enhance skin tight and elasticity, let the skin in the best state for a long time, restore the luster and vitality of the skin.
- When used, touch the back of the fuselage ion induction strip, and form a loop with the human body.
- Each mode of intelligent setting 4 minutes of service time. The instrument will enter the standby state after 4 minutes no operation and turn off automatically when 1 minute latter.

P2

SPECIFICATION ▶

Model Number:	D020	Temperature:	42°C±3°C
Vibration Frequency:	10000 rpm	Material:	ABS
Product Size:	158*55*40mm	Voltage:	3.7V
Net Weight:	88g	Power:	≤3.5W
Charging Time:	≈2H	Country of origin:	China


USER'S GUIDANCE ▶

1. Clean skin mode (Cleansing) to be used with make up ring

This mode can effectively clean the dirt in pores and make skin clear and dry



1. Take off the make up ring, attach the cotton, pour on the makeup water.

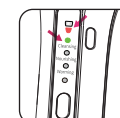

2. Long press 3s to start, clean skin brightening mold is green light. And you can choose vibration frequency strong or weak according to personal needs. vibration frequency indicator light is red. (default cleaning and brightening mode)



3. Hand contact with the back of the ion-sensing bar, according to the arrow in the picture indicating the direction of the face in a circle sliding clean.


P3

2. Nourishing and importing mode (Nourishing)

This model can effectively stimulate and relax the skin, help skin care products deeply moisturize the skin.



1. Apply daily toner or lotion to skin evenly.

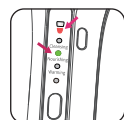

2. Long press 3s to start, select "nourish import" mode, indicator light is green. You can choose gear strength according to personal needs, vibration frequency indicator light is red.



3. Hand contact with the back of the ion-sensing bar, as shown in the arrow direction in the face from the bottom-up pull.


3. Warm heat moisturizing mode (Warming)

This model gently awakens and relaxes the skin and helps to restore its luster and vitality.



1. Apply a mask to your face after cleansing.

P4

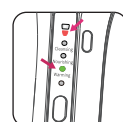

2. Long press 3s to start, choose "warm heat moisturizing" mode, indicator light is green, you can choose gear strength according to personal needs, vibration frequency indicator light is red.



3. Hand contact with the back of the ion-sensing bar, as shown in the arrow direction in the face from the bottom-up pull.

Nursing Model Description (water light muscle nurturance)

Note:

- Lead out the cleaning will the pores deep dirt to the skin surface, and then clean the face with a cleanser.
- If the finger does not contact the back of the ion induction strip, there will be no vibration response.
- This product uses the micro-current ion technology, has the slight electric flu is the normal phenomenon, please rest assured use.

Morning wake up pattern (08:00)

Arousing skin soothing edema and applying makeup

1. Clean the face
2. Apply skin care water
3. Selective nourishing mode
4. Makeup

Evening repair mode (21:00)

To nourish and nourish the living muscle floor throughout the night

1. Clean face (Cleansing mode)
2. Apply skin care (Nourishing mode)
3. Mask application (Warming mode)

P5

Machine Maintenance

- For maintenance, please check if the power is off. If the power is turned on, it is possible to use wet towels and paper towels for maintenance.
- Do not wash the fuselage, the instrument will fail, please wipe with dry cloth or soft cloth with a little soap solution.
- When there is dirt in the cotton ring, please remove it and rinse it with clean water.
- Do not wipe the fuselage with corrosive solution, such as alcohol, water, gasoline, armour removal water, etc. may cause product failure or parts cracking, discoloration, etc.

CAUTIONS ▶

1. Please read this instruction carefully before use and use this product correctly.
2. Do not clean with water. Please turn off the power before cleaning. Wipe with clean cloth or dip in a little alcohol.
3. Do not fall or strike strongly on this product, do not disassemble or modify without authorization, we do not provide free warranty service due to malfunction caused by human factors.
4. Do not approach the strong magnetic field, do not directly in the sun, do not contact with corrosive liquid, lest the machine work properly and damage.
5. Do not put this product next to high-temperature or high-power appliances to avoid fire, electric shock and other accidents.
6. This product should not stay in the same position during the use of the massage head, need to continue to move the position.
7. For the first time (first use or unused for a long time), be sure to charge, the charge indicator light from flashing to normal after the electricity is full, can be used normally.

P6

8. Non-professional do not remove this product, if the product malfunction, or not normal operation, do not continue to use, it should be sent back to the designated after-sales service intermediate maintenance.

9. Please do not allow children or limb movement restrictions, the installation of pacemakers, cardiac bypass and other people using medical equipment to use this product.

Warning! !

1. Please stop use immediately when abnormal fever, lest cause skin damage or fire.
2. Please do not fall or impact the instrument to avoid causing trouble.
3. Please do not use the power cord when it is broken, in order to avoid short-circuit and fire.

4. Please do not use the following groups:

Infectious patients, heart disease patients, high fever patients, malignant tumors, acute patients, Medical appliance, skin disease patient, pregnant/lactating person, excessive sunburner, muscle Particularly sensitive, vasodilator, Being hospitalized.

P7

GUARANTEE ▶

1. The free maintain of the guarantee card only use for the caused by product quality problems.
2. One year guarantee from the date of purchase.
3. Free to change broken parts,shipping cost is bear by customers.
4. If damage to the instrument due to force majeure is not covered by the warranty, repair service may be provided, but additional repair costs are required.
5. Please fill the gurantee clearly,without any revision. Please take good care of Gurantee card.If lost,won't be re-provided.

client's name	Contact person	
client's address	Contact Telephone	
Products name	Product model	
Purchase date	serial number	
fix records	date	Causes and handling of faults

This product has been patented, and copying is not permitted.

Patent No.: 201830385359.4

Implementation of standards: GB 4706.1-2005 / GB 4706.15-2008

P8