


Cell Body J5

USER GUIDE

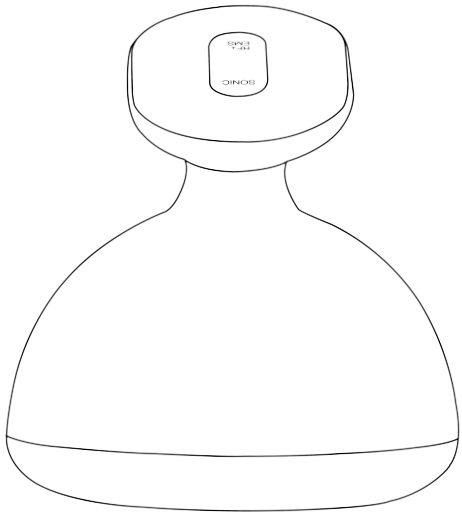

You can download the user manual in Chinese and English from the QR code.


您可以从二维码下载中英文的使用说明书。



中国語・英語の取扱説明書はQRコードよりご確認できます。

製品をご購入いただき、誠にありがとうございます。
安全を確保するため、本製品使用前に取扱説明書を十分にお読みください。
取扱説明書に保証書がついているので、本取扱説明書を大切に保管ください。

1. Name of parts

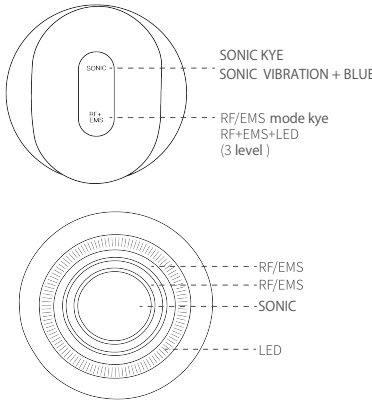

2. Product parameters

Product name	Cell Body J5
Produce Model	CX-8040
Product Size	84×81×81 mm
Battery Parameter	7.4V 1000mAh
Charging Parameter	DC5V 1A
Rated Power	8W
Ultrasonic frequency	500KHz
Material	ABS PC SUS304
Waterproof	IPX7 grade
Charging time	2H


3. Using notice

1. Please clean the part that you want to nurse before using the machine to prevent dirt from entering the skin with ultrasonic/radio-frequency or the penetration of heat energy.
2. In order to avoid the burden to body, this product is equipped with intelligent reminder (care the same area for 20 seconds and vibration stop for 1 second)to remind you to move to other areas.
3. It is recommended to use with gel or essential oil and it is a normal phenomenon that there will be a slight tingling and warm sensation while using.
4. It is recommended to try from low energy mode at first time use and then adjust it according to individual condition.
5. Three times a week is recommended. Don't nurse the same part of the body for more than 10 minutes at a time.

4. Using method


- 1.Press Sonic button or RF+EMS button to start up the device directly.
- 2.Select sonic or RF+EMS mode
sonic mode: press sonic key to start up the device instantly and turn it off by pressing again.
RF+EMS mode: Gradually strengthen the intensity by pressing the RF+EMS key there are totally 3 intensity option and turn the machine off by pressing for 4 time constantly.
The functions of Sonic+RF+EMS can be combined by pressing these function keys at the same time. The vibration can be stopped by pressing Sonic key for 3 seconds.
- 3.Please apply body gel or essential oil on the round contact.
- 4.Please stick the machine to the required slimming parts then massage your body in circles or W-curve.

5. How to use

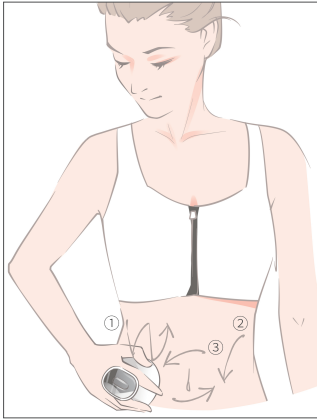

1.Waist+abdomen

Massage your waist and abdomen in W-curve way.
Massage your abdomen in circles.
Nurse for 10 minutes at a time

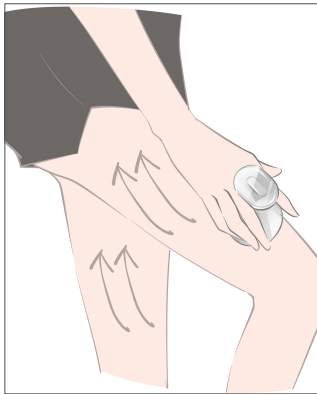


2.thighs and calves

Push the machine in flat from down to up, left legs and right legs for 5 minutes each.


3. Hip

Slide from the bottom of the thigh to the hip for 10 minutes at a time.

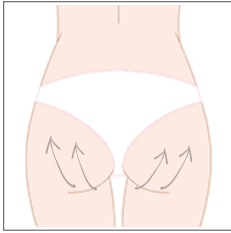


4. Arms

Massage your left and right hands for 5 minutes each from out to in.

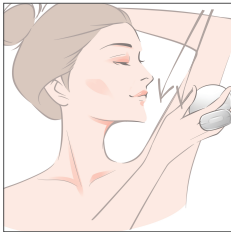


5. Chest

Lift from the back under the armpit to the chest, left and right sides for 5 minutes each.

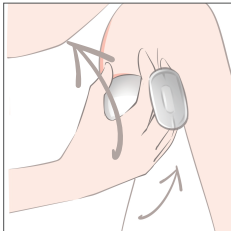

6. Neck

Pull the shoulders from down to up, left and right side for 5 minutes each.

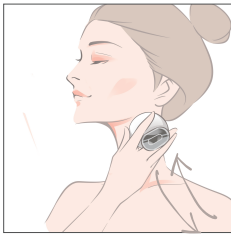

6. Mode Instruction

- 1.Red light-Skin rejuvenation mode
Press RF+EMS key to start up the skin rejuvenation mode.
Fat burning +skin tendering with red light
(EMS micro-current can activate cells and RF can accelerate the metabolism of cellulite)
It can be used on chest and body.
2. Blue light - Fat burning mode.
Press sonic key to turn on the ultrasonic fat burning mode
Fat cells burning+Blue light anti-acne(blue light is sterilizing)+Deep beauty essence lead-in.
It can be used on arm and abdomen.
3. Purple light-Fat blasting mode.
Press sonic and RF+EMS key at the same time to enter fat blasting mode.
Ultrasonic fat melting+Electric muscle stimulation(EMS) Radio-frequency(RF) micro-current muscle exercise+Vibration massage
It can be used on the thicker fat parts such as thigh/abdomen/hip etc.
Notice: this produce can't achieve the effect of breast enhancement but it can help to relieve chest problems such as enlargement and sagging.
Warm notice: please don't use it on heart area.

7. Indicator Light status

- 1.Low energy: The indicator light blinks three times then the machine will shut down.
- 2.When it is charging: the indicator light is blinking.
- 3.when it is full charged: The indicator light keeps lighting.

8. Clean and maintain

1. Please be sure to clean the gel or essence oil on the electrodes after using then store it in storage bag.
2. This product is waterproof with IPX7 grade that can be washed in water directly.
3. Please keep it in the dry environment if it is not used it for long time.

9. Notice Items

1. It is forbidden to soak for a long time.
- 2.Please confirm with your doctor if you ever have had plastic surgery. Please do not use it if you are drunk, having a cold, fever, abnormal blood pressure and other **physical discomfort, heart disease, serious disease, cancer, tuberculosis, blood disorders, allergies, dermatitis** etc.
3. Please don't use this product during the period of .pregnancy, menstruation or lactation.
- 4.Please don't use it on skin area of inflammation and around eyes, lips, vocal cords and other sensitive areas and don't use it on the same part for long time.
5. This product is for beauty care, please don't use in other ways and don't disassemble or change or repair it by yourself.
6. Please inspect the device before using. If the contact is cracked, please stop using it.
- 7.Please stop using if you feel uncomfortable.

10. Common Issues

Common issues	Solutions
If it is unable to work after starting up	Check if the power is low and then restart it after fully charged
If it is unable to charge	1.Check if the charging cable is well connected. 2.Check if the charging port is blocked. 3.Please contact customer service and send the product to repair if it is still unable to work.
Waterproof	This product is waterproof with IPX7 grade which can be used for nursing while taking in a bath or shower but please don't use it with bath liquid or taking hot spring bath or sauna.

11. Warranty instruction

- 1.Please keep this warranty card carefully as maintenance certificate.
- 2.The warranty period is 12 months from the date of purchase.
- 3.We will offer the fixing service and replace components for free once we checked out
- 4.product is abnormal within **the warranty period under the condition of** normal using and using and maintenance

The followings are not covered by the warranty

- 1.Product is broken due to improper operation.
- 2.**User has fixed,changed,repared or replaced any component of the product.**
- 3.The serial number of machine has been altered.
- 4.Damage caused by careless use or leakage of water or other substances into the machine.