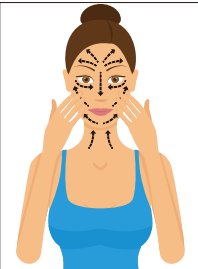
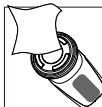
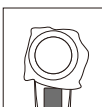
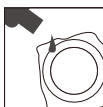
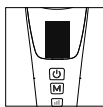


















<div>Skin rejuvenation instrument</div> <div>SKIN BEAUTY UP</div> <div>User Manual▶</div>	<div>Product structure</div> <div><div><div><div><div>Alloy introduction head</div><div>LED lights</div><div>Ring for fix cotton</div></div><div><div>Display</div><div>Power switch</div><div>Mode adjustment</div><div>Level</div></div><div><div>Ion induction strip</div><div>USB charging</div></div></div></div></div> <div>-1-</div>	<div>Specification</div> <table><tr><td>Name</td><td>SKIN BEAUTY UP</td></tr><tr><td>Model</td><td></td></tr><tr><td>Charging voltage</td><td>DC 5V</td></tr><tr><td>Rated power</td><td><5W</td></tr><tr><td>Rated voltage</td><td>About 3.7V</td></tr><tr><td>Charging time</td><td>About 2.5H</td></tr><tr><td>Vibration frequency</td><td>10000(±10%) RPM</td></tr><tr><td>Temperature range</td><td>42±5℃</td></tr><tr><td>Product size</td><td>175*55*67mm</td></tr></table> <div>Machine maintenance</div> <div><div>● During maintenance, please confirm whether the power is turned off. If the power is turned on, or if maintenance is performed with wet towels and tissues, there may be electrical stimulation.</div><div>● Do not wash the body with water, the instrument will malfunction, please wipe with a dry cloth or a soft cloth moistened with soapy liquid.</div><div>● When the cotton pad retaining ring is dirty, remove it and rinse it with clean water.</div><div>● Do not wipe the body with corrosive solutions, such as Organic solvent, gasoline, nail polish remover, etc., may cause product failure or cracking and discoloration of accessories.</div></div> <div>-2-</div>	Name	SKIN BEAUTY UP	Model		Charging voltage	DC 5V	Rated power	<5W	Rated voltage	About 3.7V	Charging time	About 2.5H	Vibration frequency	10000(±10%) RPM	Temperature range	42±5℃	Product size	175*55*67mm	<div>Recommended massage method</div> <div><div><div></div><div><div>Chin: Place the massage head in the middle of the chin, and use the arc shape to Chin lifting at an average speed of 2s per time.</div><div>Interfacial line: Place the fuselage at the lower side of the temples, and massage to the jawbone at an average speed of 2s per time from bottom to top.</div><div>Cheeks: Place the massage head next to the corner of the lips, in the direction of the earlobe, in a curved manner, and massage at an average speed of 2s per time.</div><div>Eyes: Massage around the eye circle at average speed of 2s per time.</div></div></div><div><div>8:00</div><div>Morning Rejuvenation Mode</div><div>Awaken the skin + soothe edema + apply makeup</div><div>Clean face</div><div>Apply toner</div><div>Make up</div><div>21:00</div><div>Evening repair mode</div><div>Revitalize the bottom of the muscle + strengthen nutrient absorption + overnight repair</div><div>Clean face</div><div>Apply toner</div></div><div>-3-</div></div>	<div>The mode and function of LED colorful light</div> <div><div>Please switch on to the second mode photo rejuvenation, and then press the level key, according to personal skin needs with six kinds of colored light to care for the skin. You can use LED mode together with skin care product.</div><div>Red light: promote blood circulation, anti-wrinkle and anti-aging (wavelength 700 ± 10nm)</div><div>Blue light: calm the skin and inhibit the growth of acne (wavelength 475 ± 10nm)</div><div>Yellow light: Regulate cell activity and diminish freckles (wavelength 600 ± 10nm)</div><div>Green light: remove dark spots and yellow spots, whiten (wavelength 550 ± 10nm)</div><div>Orange light: diminish erythema and redness of pigment (wavelength 625 ± 10nm)</div><div>Purple light: calms and reduces inflammation, suitable for oily skin (wavelength 420 ± 10nm)</div></div> <div>-4-</div>
Name	SKIN BEAUTY UP																					
Model																						
Charging voltage	DC 5V																					
Rated power	<5W																					
Rated voltage	About 3.7V																					
Charging time	About 2.5H																					
Vibration frequency	10000(±10%) RPM																					
Temperature range	42±5℃																					
Product size	175*55*67mm																					
<div>Instructions</div> <div><div>How to use clean mode</div><div>This mode can effectively clean the dirt in the pores and remove makeup</div><div><div><div></div><div></div><div></div></div><div><div></div><div></div></div></div><div>The above steps are used to remove makeup. If you did not apply makeup, you can replace the makeup remover with toner to clean your skin.</div><div><div>How to use photo rejuvenation mode</div><div>This mode is to introduce nutrition essence into the dermis layer to make the skin more young and moisturize</div><div><div><div></div><div></div></div><div><div></div><div></div></div></div></div><div>-5-</div></div>	<div><div>How to use firming and lifting mode</div><div>This mode can enhance cell activity, lanti-warinkle and help skin products to nourish deeply</div><div><div><div></div><div></div></div><div><div></div><div></div></div></div><div><div>How to use eye care mode</div><div>This mode relieves eye fatigue and helps lighten dark circles and anti crows feet</div><div><div><div></div><div></div></div><div><div></div></div></div></div><div>-6-</div></div>	<div>Function</div> <div><div><div></div><div><div>CLEAN DETOXIFICATION</div><div>High-frequency vibration + heating to open the pores, heating to open the pores, shattering the dirt and exporting the deep-seated garbage at the deep of pores, clearing the pores and reducing blackhead formation.</div></div></div><div><div></div><div><div>PHOTOREJUVENATION</div><div>Fully introduce skin care products, warm massage to promote blood circulation, improve skin firmness and elasticity, and keep skin in the best condition for a long time. Color light care enhances cell activity, achieves whitening and skin rejuvenation, and anti wrinkles. Comprehensive restorative skin radiance and vitality.</div></div></div><div><div></div><div><div>FIRMING AND LIFTING</div><div>The micro current can stimulate facial muscles through skin tissues by activating gentle electric waves, activate cell regeneration, restore skin elasticity, let the face get exercise, and also repair the dermis layer to achieve the effect of firming the skin.</div></div></div><div><div></div><div><div>EYE CARE</div><div>Warm and sonic vibration assists eye care can soothes eye fatigue and works better with eye cream. Compared with manual surface application, it can more effectively introduce eye cream into the deep cell, leaving the skin around the eyes more moisture and young.</div></div></div></div> <div>-7-</div>	<div>Safety Precautions</div> <div><div>1. Please read the instructions carefully before use, and then use the product correctly in accordance with the operation sequence.</div><div>2. Please do not disassemble or modify the product by yourself. Please consult a special repair shop if the product fails.</div><div>3. This product is only used for beauty care aids, please do not use it for purposes other than instructions.</div><div>4. Please don't use it under high temperature and high humidity environment, so as not to lose or damage the function of the instrument!</div><div>5. Do not place this product near high-temperature or high-power electrical appliances to avoid fire and electric shock accidents.</div><div>6. Please do not charge in a humid environment to avoid short circuit, electric shock, fire, etc. of the instrument.</div><div>7. Please do not let children or people with limited physical activity use, minors need to use under adult supervision.</div><div>8. Massage head can be wiped clean with a cleaning cloth or wiped with a little alcohol for disinfection.</div><div>9. This product has no waterproof function, please do not place it in water or in an environment with corrosive chemicals.</div><div>Caveat</div><div>1. Please stop using it immediately when there is abnormal heat to avoid scalding and fire.</div><div>2. Please don't drop or hit the instrument to avoid instrument failure.</div><div>3. The essence needs to be supplemented during use to avoid dry skin.</div><div>4. When the power cord is damaged, please do not use it to avoid short circuit and fire.</div><div>5. Please do not use the following people: People with heart disease, people with high fever, people with infectious diseases, people with malignant tumors, people with skin wounds, acute diseases. People who are using medical appliances, patients with skin diseases, pregnant women who are breastfeeding, people with excessive sunburn, people with particularly sensitive skin, people with vasodilators, and those who are in hospital.</div></div> <div>-8-</div>	<div>Warranty card</div> <div><div><div>保証書</div><div>品名SKIN BEAUTY UP</div><div>お名前</div><div>ご住所 〒</div><div>TEL</div><div>ご購入日 年 月 日 保証期間/ご購入日より1年間</div><div>販売店印</div></div><div><div>発売元 株式会社 エヴァンジェール 製造販売元 株式会社 アグレックス</div><div>■本社 千321-0164 ■本社 千361-0023</div><div>栃木県宇都宮市双葉 1-13-50 千葉県千葉市美浜区中瀬 1-3 MTG B 棟 9F</div><div>TEL 0120-169-119 TEL 043-350-0111 FAX 043-350-0114</div></div></div> <div>-9-</div>																		